


# Evetric®

1

Self-Etch

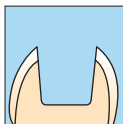


or

Selective Enamel-Etch






Total Etch

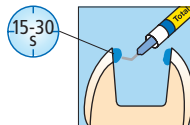

2

Evetric® Bond

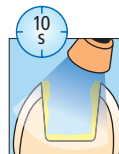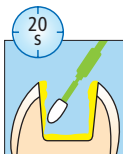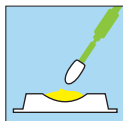

3

Optional  
Tetric EvoFlow®





Bluephase® Style

$\geq 500 \text{ mW/cm}^2$	20 s
$\geq 1'000 \text{ mW/cm}^2$	10 s
Bluephase® Style	10 s

4

Evetric®

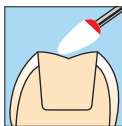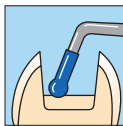

max. 2 mm  
Dentin: max. 1.5 mm

CE 0123



See Instructions

Date information prepared: 2015-04-24  
673606/Rev. 0

ivoclar  
vivadent  
clinical

A1



A2



A3



A3.5



B2



B3



C2



T



A3.5  
Dentin



B2  
Dentin

