

Compoglass® F

Compomer restorative material



A
combination
of
glass ionomer and
composite



ivoclar
vivadent
passion vision innovation

Compoglass® F

Compomer restorative material




Compoglass® F is a light-curing compomer-based restorative material.

The product has been successfully used by dentists for more than 15 years.

Compoglass F combines the benefits of glass ionomers with those of light-curing composite resins. The product contains aluminium fluorosilicate glass and ytterbium trifluoride and therefore releases high amounts of fluoride. Due to advanced composite technology, the material exhibits high strength and wear resistance.

Compoglass F is used to fill cavities in deciduous teeth and in cervical areas.

High fluoride release

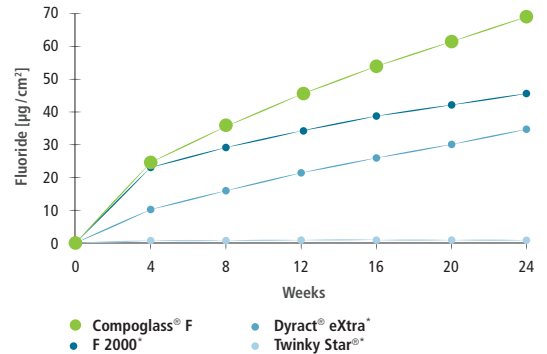

Source: R&D Ivoclar Vivadent AG, Schaan, 2014
* These brands are not registered trademarks of Ivoclar Vivadent AG.

Compoglass® Flow is the flowable version of Compoglass F.

The ultra-thin and long Cavifil tips allow the material to be syringed into the cavity quickly and easily. As a result of its excellent flow properties, the material readily adapts to the cavity margins.

The universal adhesive Adhese® Universal can be used in combination with these materials to streamline the treatment process. Adhese Universal is applied directly in the mouth from the unique and convenient VivaPen® delivery form.



Delivery forms

Compoglass F

574992AN	Intro Pack		
546980AN	A1	Refill 20 x 0.25 g	546983AN A3.5 Refill 20 x 0.25 g
546981AN	A2	Refill 20 x 0.25 g	546985AN A4 Refill 20 x 0.25 g
546982AN	A3	Refill 20 x 0.25 g	546979AN 105 Refill 20 x 0.25 g

Compoglass Flow

550167AN	A2	Refill 20 x 0.25 g
550168AN	A3	Refill 20 x 0.25 g



This product forms a part of the "Direct Restoratives" product category. Products from this category are optimally coordinated with each other.

Ivoclar Vivadent AG
Benderstr. 2
9494 Schaan
Liechtenstein
Tel. +423 235 35 35
Fax +423 235 33 60
www.ivoclarvivadent.com

ivoclar
vivadent®
passion vision innovation